

# PEDI-PAK<sup>®</sup>

## Cures Your Splitting Headache!

**NOW!**

- Quad or Single
- Needleless Y Site<sup>™</sup> or Standard Injection Site
- Latex-free

PEDI-PAK<sup>®</sup> 75 mL Transfer Packs are the ideal solution for splitting an adult size unit of blood into four equal pediatric size units. And, for maximum patient safety, each PEDI-PAK<sup>®</sup> has a completely sterile, pyrogen free fluid path.

- Maximum usage of adult size blood units
- Reduces risk of antigen exposure to neonates
- Ideal for preparing small units of CMV negative and rare blood or FFP
- Compatible with sterile connecting technology
- Meets the small volume storage requirements for the blood bank or lab
- Latex-free

Distributed Exclusively By

**NATIONAL HOSPITAL SPECIALTIES**

*Blood Products Specialties LLC division*

65 Commerce Way, Hackensack, NJ 07601

Our Needleless "Y" Site<sup>™</sup> complies with OSHA's "Federal Needlestick Safety and Prevention Law."

**PEDI-PAK<sup>®</sup> SINGLE** ▶



**CALL: 800-344-6087/201-488-1174 FAX: 201-488-0983 WEBSITE: [www.nhsusa.com](http://www.nhsusa.com)**

# PEDI-PAK<sup>®</sup> 75 mL Transfer Packs

## NEW FEATURES


- NEW!** ■ Available with Needleless "Y" Site<sup>™</sup> or Standard Injection Site
- NEW!** ■ Now latex-free
- Quad or Single

## SPECIAL FEATURES


- Four interconnected transfer packs or single transfer bag (6.4 cm x 13.3 cm) with a 75 mL capacity
- Convenient access port on each aliquot
- With or without an in-line filter (195 microns)
- Completely sterile, pyrogen free fluid path
- Compatible with sterile connecting technology
- Fully protected spike for safe handling
- Sterile packaging for convenient storage
- Easy-to-follow instructions
- Needleless "Y" Site<sup>™</sup> helps avoid needle sticks

## SUGGESTED APPLICATIONS


- For preparing whole blood/red cell pediatric aliquots
- For preparing aliquots of single donor fresh frozen plasma (FFP)
- For preparing small volume units of rare and CMV negative blood
- For storage of small volumes of research specimens and other special products
- For use with a pediatric syringe filter

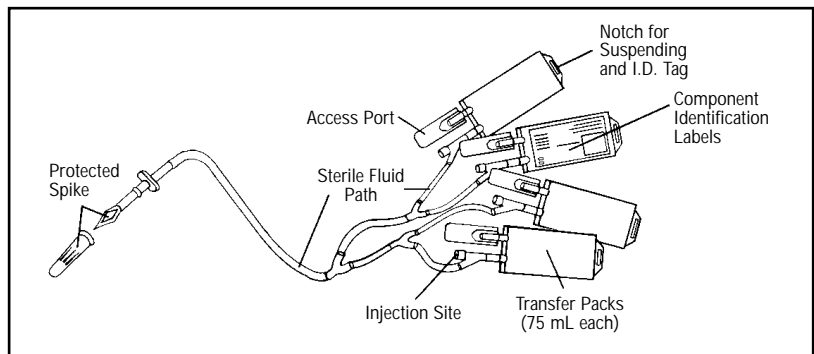

PEDI-PAK<sup>®</sup> QUAD 75 mL Transfer Pack  
Cat # 402-04      Packaged: 36 per case

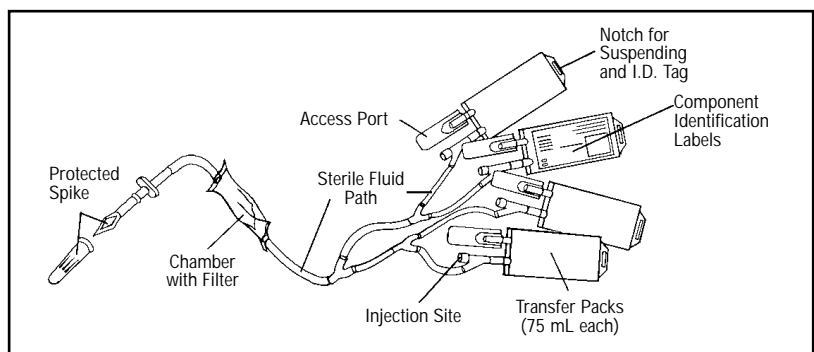

PEDI-PAK<sup>®</sup> QUAD 75 mL Transfer Pack with In-line Filter  
Cat # 402-04F      Packaged: 36 per case

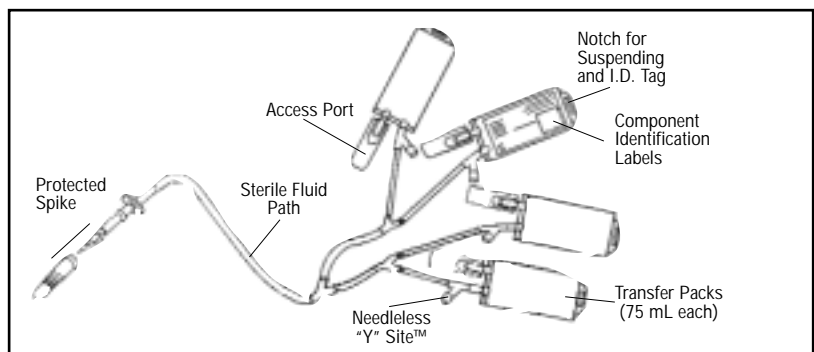

PEDI-PAK<sup>®</sup> QUAD 75 mL Transfer Pack with Needleless "Y" Site<sup>™</sup>  
Cat # 402-024      Packaged: 36 per case

